

Manuals+

User Manuals Simplified.



Berrcom JXB-178 Non-Contact Infrared Thermometer User Guide

[Home](#) » [Berrcom](#) » Berrcom JXB-178 Non-Contact Infrared Thermometer User Guide 




Links Medical Products Inc.[®] JXB-178
Non-Contact Infrared Thermometer
Quick Start Guide

Contents [[hide](#)]

- [1 USING THE THERMOMETER](#)
- [2 IMPORTANT: BEFORE TAKING THE FIRST MEASUREMENT](#)
- [3 CHANGING TEMPERATURE MODE](#)
- [4 CHANGING TEMPERATURE UNIT](#)
- [5 TAKING A MEASUREMENT](#)
- [6 LOW BATTERY DETECTION](#)
- [7 CHANGING THE BATTERIES](#)
- [8 TROUBLESHOOTING](#)
- [9 File Downloads](#)
- [10 Related Manuals](#)

USING THE THERMOMETER

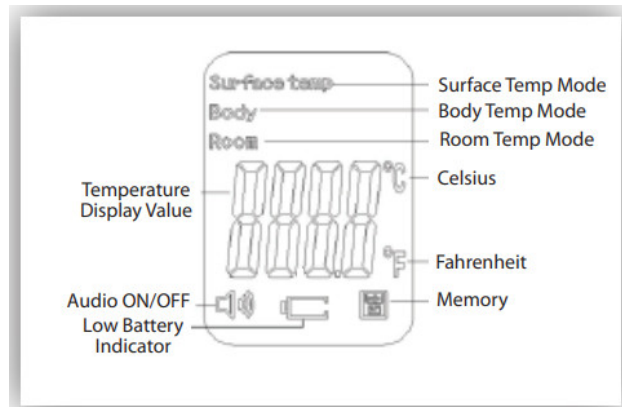

IMPORTANT: BEFORE TAKING THE FIRST MEASUREMENT

The thermometer needs to acclimate to the ambient environment (temperature) in which measurements will occur. It is important to allow the JXB-178 to acclimate to the room temperature for at least 15 to 20 minutes prior to first using the thermometer.

Subsequent usage – each time there is a significant change in the ambient temperature due to a change in the environment, allow the JXB-178 to acclimate to this ambient temperature for at least 15 minutes prior to using it.

CHANGING TEMPERATURE MODE

While the thermometer is turned on, press the MODE button to toggle through the three temperature mode options: BODY, SURFACE or ROOM.

NOTE: When first using the thermometer, the default is set to BODY mode.

NOTE: Measurement ranges for each temperature mode are as follows:

BODY mode: 89.6°F to 109°F (32°C to 42.9°C)

SURFACE mode: 32°F to 140°F (0°C to 60°C)

ROOM mode: 50°F to 104°F (10°C to 40°C)

NOTE: Once the thermometer is powered on, the last temperature mode measured will default as the current temperature mode.


NOTE: Always use the BODY mode to measure internal human temperature.

CHANGING TEMPERATURE UNIT

While the thermometer is turned on, press the MODE button for 2 seconds, the screen will display “F1”, then press the MODE button to transfer between degrees Celsius and Fahrenheit, confirm by pressing the MEM button.


TAKING A MEASUREMENT

1. Press the ON/SCAN button to power on the thermometer.
When the device is ready, the LCD screen will turn green.
2. Position the thermometer between 1.2” to 2” (3 to 5 cm) from the patient’s forehead.
3. Press the ON/SCAN button to take measurements.
4. Measurement will be displayed on the LCD screen.
 - To repeat the process, return to step 2.
 - Allow 3 to 5 seconds between measurements.
5. The thermometer automatically turns off after 30 seconds of inactivity.

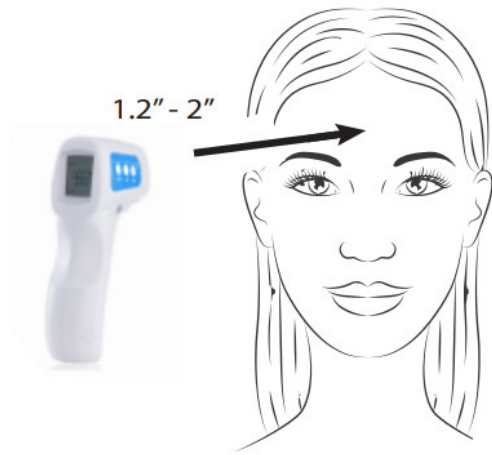


Note: High-temperature alarm default setting is 100.4°F (38°C).

LOW BATTERY DETECTION

If the battery voltage is not more than 2.7V, the JXB-178 Non-contact Infrared Thermometer will display the flashed symbol to  remind you to replace both batteries.

CHANGING THE BATTERIES

The JXB-178 Non-contact Infrared Thermometer requires 2 AA alkaline batteries. When the LCD screen displays , both batteries should be replaced. Open the lid on the bottom of the thermometer and replace the batteries, noting the correct battery positioning. Wait 15-20 minutes for the unit to acclimate to the ambient environment (temperature) before performing a test.

TROUBLESHOOTING

Reasons for Lo message display

Advice

The temperature reading is hampered by hair, perspiration or makeup.

Make sure there is no obstruction or dampness prior to taking temperature.

Temperature hampered by an air draft or dramatic change in ambient temperature.

Make sure there is no air blowing in the area of use, this could affect the infrared reading.

Temperature readings are too close together and the thermometer did not have time to reboot.

Pause for 3-5 seconds minimum between two readings; a 5-second pause is recommended.

Measuring too far from the forehead. The ambient temp < 30°F (0°C).

Take measurements at the recommended distance of 1.2" to 2" (3 to 5 cm) from the forehead.

Measuring after cold compress, taking medicine or bathing.

Make sure the measuring environment temperature is higher than 32°F (0°C). If any of these apply, make sure to wait 30 minutes and test again.

Reasons for Hi message display

When body temperature is above 109.2°F (42.9°C).

F2 value is set higher than desired.

The room/ambient temperature is above 104°F (40°C).

The surface temperature exceeds 140°F (60°C).

Advice

Temperature reading exceeds measurement range of 89.6°F to 109.2°F (32°C to 42.9°C). Make sure there is no hot air blowing in the area of use. If necessary, immediately seek medical attention.

Reset high alarm temperature setting closer to or at the default value of 100.4°F (38°C).

Make sure to operate thermometer between 50°F to 104°F (10°C to 40°C).

Temperature reading exceeds measurement range of 32°F to 140°F (0°C to 60°C).



LMP Item# LMP010



1 Year Limited Warranty **Links Medical Products Inc.®**

9247 Research Drive
Irvine, CA 92618 USA
888.425.1149
www.linksmed.com

File Downloads



Berrcom User Guide [Download](#) [optimized]
JXB-178, Non-Contact Infrared Thermometer, JXB-178 Non-Contact Infrared Thermometer [Download](#)

Related Manuals

1. [Berrcom Non-Contact Infrared Thermometer JXB-178 User Guide](#) Non-contact Infrared Thermometer Model:JXB-178 Quick guide Measure Temperature Indoors Before...
2. [Berrcom Non-Contact Infrared Digital Thermometer JXB-178 User Manual](#) Berrcom Non-Contact Infrared Digital Thermometer JXB-178 THE MANUFACTURER RESERVES...
3. [Berrcom JXB-183 Non-Contact Infrared Thermometer User Guide](#) JXB-183 Non-Contact Infrared Thermometer Quick Start Guide USING THE THERMOMETER...
4. [Berrcom JXB-178 Non-Contact Infrared Thermometer User Manual](#) JXB-178 Non-Contact Infrared Thermometer...
5. [VERIDIAN JXB-311 Non-Contact Infrared Thermometer User Guide](#) Non-Contact Infrared Thermometer Model JXB-311 / Ref No 09-313 These...
6. [Berrcom Infrared Thermometer User Guide](#) Berrcom JXB-178 Non-Contact Infrared Thermometer Quick Start Guide USING THE...